













Commonly Used Symbols











Glossary of symbols commonly used on
DeRoyal Medical Device Packaging and Labeling

symbol	explanatory text	standard reference	standard title
	Indicates the manufacturer's catalogue number so that the medical device can be identified.	ISO 15223-1, Clause 5.1.6	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Catalogue or model number		ISO 7000- 2493	Graphical symbols for use on equipment.
	Indicates the manufacturer's batch code so that the batch or lot can be identified.	ISO 15223-1, Clause 5.1.5	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Batch code		ISO 7000-2492	Graphical symbols for use on equipment.
	Indicates the manufacturer's serial number so that a specific medical device can be identified.	ISO 15223-1, Clause 5.1.7	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Serial number		ISO 7000-2498	Graphical symbols for use on equipment.
	Indicates a medical device that needs to be protected from moisture.	ISO 15223-1, Clause 5.3.4	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Keep dry		ISO 7000-0626	Graphical symbols for use on equipment.
	Indicates a medical device that needs protection from light sources.	ISO 15223-1, Clause 53.2	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Keep away from sunlight		ISO 7000-0624	Graphical symbols for use on equipment.
	Indicates the medical device manufacturer.	ISO 15223-1, Clause 5.1.1	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Manufacturer		ISO 7000-3082	Graphical symbols for use on equipment.
	Indicates the date when the medical device was manufactured.	ISO 15223-1, Clause 5.1.3	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Date of manufacture		ISO 7000-2497	Graphical symbols for use on equipment.
Rx Only	Requires prescription in the United States.	21 CFR 801.15(c)(1)(i)F	Labeling-Medical devices; prominence of required label statements.
		21 CFR 801.109	Labeling-Prescription devices.
Prescription only			
	Indicates a medical device that should not be used if the package has been damaged or opened.	ISO 15223-1, Clause 5.2.8	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
		ISO 7000-2606	Graphical symbols for use on equipment.
Do not use if package is damaged			
	Indicates the presence of natural rubber or dry natural rubber latex as a material of construction within the medical device or the packaging of a medical device.	ISO 15223-1, Clause 5.4.5	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Contains natural rubber latex			
	Indicates the materials of construction of the device or its packaging were not manufactured with natural rubber or dry natural rubber latex.	ISO 15223-1, Clause 5.4.5, Reference Annex B for negation symbol	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Not made with natural rubber latex			



Commonly Used Symbols

Glossary of symbols commonly used on
DeRoyal Medical Device Packaging and Labeling

symbol	explanatory text	standard reference	standard title
	Indicates the date after which the medical device is not to be used.	ISO 15223-1, Clause 5.1.4	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
Use by		ISO 7000-2607	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
	Indicates the need for the user to consult the instructions for use.	ISO 15223-1, Clause 5.4.3	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
Consult instructions for use		IEC 60601-1, Table D.1, Symbol 11	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
		ISO 7000-1641	Medical electrical equipment — Part 1: General requirements for basic safety and essential performance.
	Indicates the need for the user to consult the instructions for use for important cautionary information such as warnings and precautions that cannot, for a variety of reasons, be presented on the medical device itself.	ISO 15223-1, Clause 5.4.4	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
Caution		IEC 60601-1, Table D.1, Symbol 10	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
	Indicates a medical device that is intended for one use or for use on a single patient during a single procedure.	ISO 15223-1, Clause 5.4.2	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
Do not reuse		ISO 7000-1051	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
	Refer to instruction manual/ booklet.	IEC 60601-1, Table D.2, Symbol 10	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
Follow instructions for use			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
	Signifies European technical conformity.	765/2008/EC 768/2008/EC MDD 9 3/42/EEC Articles 4,11,12,17, Annex II)	Medical electrical equipment — Part 1: General requirements for basic safety and essential performance.
CE marking			The requirements for accreditation and market surveillance relating to the marketing of products; Medical Device Directive.
	Indicates the authorized representative in the European Community.	ISO 15223-1, Clause 5.1.2	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
Authorized European representative			Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
	Indicates the temperature limits to which the medical device can be safely exposed.	ISO 15223-1, Clause 5.3.7	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
Storage temperature range		ISO 7000-0632	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
	Indicates the lower limit of temperature to which the medical device can be safely exposed.	ISO 15223-1, Clause 5.3.5	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
Lower limit of storage temperature		ISO 7000-0534	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
	Indicates the upper limit of temperature to which the medical device can be safely exposed.	ISO 15223-1, Clause 5.3.6	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
Upper limit of storage temperature		ISO 7000-0533	Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.



Commonly Used Symbols


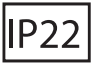







Glossary of symbols commonly used on
DeRoyal Medical Device Packaging and Labeling

symbol	explanatory text	standard reference	standard title
	Indicates the range of humidity to which the medical device can be safely exposed.	ISO 15223-1, Clause 5.3.8	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements.
Humidity Limitation		ISO 7000-2620	Graphical symbols for use on equipment.
	Indicates a medical device that is not to be re-sterilized.	ISO 15223-1, Clause 5.2.6	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements. Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Do not re-sterilize		ISO 7000-2608	Graphical symbols for use on equipment.
	Indicates a medical device that has not been subjected to a sterilization process.	ISO 15223-1, Clause 5.2.7	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements. Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Non sterile		ISO 7000-2609	Graphical symbols for use on equipment.
	Indicates a medical device that has been subjected to a sterilization process.	ISO 15223-1, Clause 5.2.1	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements. Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Product subjected to a sterilization process		ISO 7000-2499	Graphical symbols for use on equipment.
	Indicates a medical device that has been sterilized using irradiation.	ISO 15223-1, Clause 5.2.4	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements. Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Sterilized using irradiation		ISO 7000-2502	Graphical symbols for use on equipment.
	Indicates a medical device that has been sterilized using aseptic techniques.	ISO 15223-1, Clause 5.2.2	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements. Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Sterilized using aseptic techniques		ISO 7000-2500	Graphical symbols for use on equipment.
	Indicates a medical device that has been sterilized using ethylene oxide.	ISO 15223-1, Clause 5.2.3	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements. Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Sterilized by ethylene oxide treatment		ISO 7000-2501	Graphical symbols for use on equipment.
	Indicates a medical device that has been sterilized using steam or dry heat.	ISO 15223-1, Clause 5.2.5	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements. Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Sterilized using steam or dry heat		ISO 7000-2503	Graphical symbols for use on equipment.
	Indicates there are potential biological risks associated with the medical device.	ISO 15223 5.4.1	Medical Devices — Symbols to be used with information to be supplied by the manufacturer — Part 1: General Requirements. Medical devices — Symbols to be used with medical device labels, labeling and information to be supplied.
Biological Risks		ISO 7000-0659	Graphical symbols for use on equipment.
	Double insulated appliance	IEC 60417-5172	Graphical symbols for use on equipment.
Double insulated appliance			



Commonly Used Symbols

Glossary of symbols commonly used on
DeRoyal Medical Device Packaging and Labeling

symbol	explanatory text	standard reference	standard title
	Type BF applied part	IEC 60417-5333	Graphical symbols for use on equipment.
Type BF applied part			
	Classifies and rates the degree of protection provided by casing/ enclosure of electrical device.	IEC 60529*	Graphical symbols for use on equipment.
IP code, International Protection Marking	<ul style="list-style-type: none"> • Protection against solid objects > 12.5 mm; • Protection against water sprayed up to 15 from vertical 	*(does not have a reference number per 60601-1 [2005] p.645)	
	Magnetic resonance safe	ASTM F2503	Standard Practice for Marking Medical Devices and Other Items for Safety in the Magnetic Resonance Environment.
MR Safe			
	Magnetic resonance conditional	ASTM F2503	Standard Practice for Marking Medical Devices and Other Items for Safety in the Magnetic Resonance Environment.
MR Conditional			
	Magnetic resonance unsafe	ASTM F2503	Standard Practice for Marking Medical Devices and Other Items for Safety in the Magnetic Resonance Environment.
MR Unsafe			
	Waste of electrical and electrical equipment	EN 50419	Standard Practice for Marking Medical Devices and Other Items for Safety in the Magnetic Resonance Environment.
WEEE			
	Non-ionizing electromagnetic radiation	IEC 60417 REF 5140	Standard Practice for Marking Medical Devices and Other Items for Safety in the Magnetic Resonance Environment.
Non-ionizing Radiation			
	Indicates the materials of construction do not include the phthalate plasticizer DEHP	BS EN 15986 Annex B	Symbol for use in the labeling of medical devices – Requirements for labeling of medical devices containing phthalates.
Does not contain DEHP			
	Indicates the item is a medical device.	ISO 15223 – 1, Clause 5.7.7	Medical Devices – Symbols to be used with information to be supplied by the manufacturer – Part 1: General Requirements.
Medical Device			